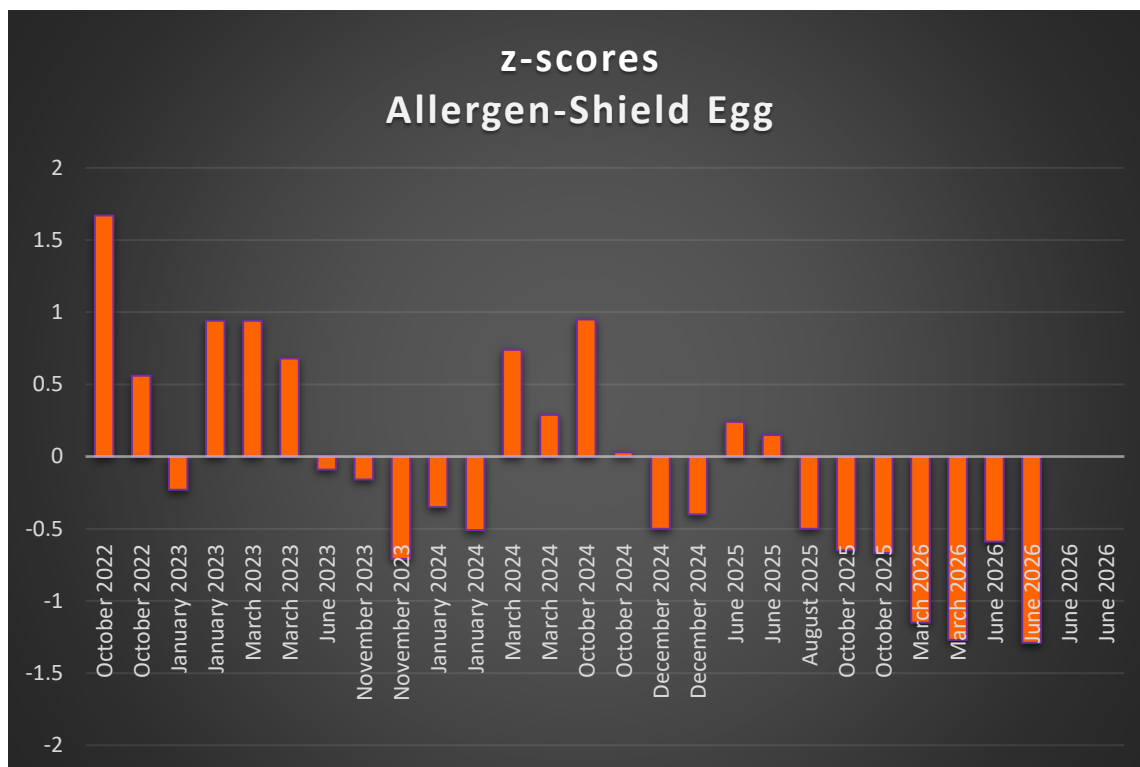




| Proficiency Test                                                                                   | Assigned value (mg/kg)   | Result (mg/kg) | Z-score |
|----------------------------------------------------------------------------------------------------|--------------------------|----------------|---------|
| FAPAS Egg in Cake Mix 27470AQC April – June 2026                                                   | 5.81(allergenic protein) | 5.82           | 0.0     |
| FAPAS Egg in Cake Mix 27470AQC April – June 2026                                                   | 12.1                     | 12.1           | 0.0     |
| FAPAS Egg in Cake Mix 27470BQC April – June 2026                                                   | NOT DETECTED             | <1.0           | -       |
| BIPEA Allergens in Industrial Process Water 77D - Industrial Process Water – 21 - 1077 – June 2026 | 10.6                     | 8.9            | -0.59   |
| BIPEA Allergens in Industrial Process Water 77D - Industrial Process Water – 21 - 1177 – June 2026 | 18.5                     | 16.3           | -1.29   |

|                                                                                                    |              |      |       |
|----------------------------------------------------------------------------------------------------|--------------|------|-------|
| BIPEA Allergens in Industrial Process Water 77D - Industrial Process Water – 21 - 1277 – June 2026 | NOT DETECTED | <0.4 | -     |
| FAPAS Gluten, Milk & Egg in Cooked Biscuit 27459 January – March 2026                              | 12.7         | 11.7 | 0.1   |
| BIPEA Food Allergens – EGG 77B – Pasta – 11 - 2577 - March 2026                                    | 42.5         | 23.7 | -1.15 |
| BIPEA Food Allergens – EGG 77B – Pasta – 10 - 2677 - March 2026                                    | NOT DETECTED | <0.4 | -     |
| BIPEA Food Allergens – EGG 77B – Pasta – 10 - 2777 - March 2026                                    | 74.9         | 58.6 | -1.27 |
| FAPAS Gluten, Milk & Egg in Cake Mix 27438B July – August 2025                                     | 18.8         | 16.3 | -0.5  |
| BIPEA Food Allergens – EGG 77B – Cake Mix – 11 - 0477 - October 2025                               | NOT DETECTED | <0.4 | -     |
| BIPEA Food Allergens – EGG 77B – Cake Mix – 11 - 0577 - October 2025                               | 24.8         | 20.5 | -0.67 |
| BIPEA Food Allergens – EGG 77B – Cake Mix – 11 - 0677 - October 2025                               | 45.4         | 40.5 | -0.65 |
| BIPEA Allergens in Industrial Process Water 77D - Industrial Process Water – 18 - 1077 – June 2025 | 66.0         | 70.4 | 0.24  |
| BIPEA Allergens in Industrial Process Water 77D - Industrial Process Water – 19 - 1177 – June 2025 | 38.4         | 39.9 | 0.15  |
| BIPEA Allergens in Industrial Process Water 77D - Industrial Process Water – 19 - 1277 – June 2025 | NOT DETECTED | <0.4 | -     |
| FAPAS Gluten, Egg and Milk in Cake Mix 27400A July – September 2024                                | 11.3         | 9.87 | -0.5  |
| FAPAS Gluten, Egg and Milk in Cake Mix 27400B July – September 2024                                | 22.2         | 20.1 | -0.4  |

|                                                                                                       |              |      |       |
|-------------------------------------------------------------------------------------------------------|--------------|------|-------|
| BIPEA Food Allergens – EGG 77B – Cake Mix – 10 - 0677 - October 2024                                  | 23.5         | 23.7 | 0.03  |
| BIPEA Food Allergens – EGG 77B – Cake Mix – 10 - 0577 - October 2024                                  | NOT DETECTED | <0.4 | -     |
| BIPEA Food Allergens – EGG 77B – Cake Mix – 10 - 0477 - October 2024                                  | 32.9         | 41.0 | 0.95  |
| BIPEA Food Allergens – EGG 77B - Pasta – 09 - 2577 - March 2024                                       | NOT DETECTED | <0.4 | -     |
| BIPEA Food Allergens – EGG 77B - Pasta – 08 - 2677 - March 2024                                       | 22.2         | 24.2 | 0.29  |
| BIPEA Food Allergens – EGG 77B - Pasta – 08 - 2777 - March 2024                                       | 54.1         | 63.6 | 0.74  |
| BIPEA Allergens in Industrial Process Water 77D - Industrial Process Water – 16 - 1077 – January 2024 | NOT DETECTED | <0.4 | -     |
| BIPEA Allergens in Industrial Process Water 77D - Industrial Process Water – 16 - 1177 – January 2024 | 12.8         | 10.2 | -0.51 |
| BIPEA Allergens in Industrial Process Water 77D - Industrial Process Water – 16 - 1277 – January 2024 | 17.0         | 14.3 | -0.35 |
| BIPEA Food Allergens – EGG 77B - Cake Mix – 09 – 0477 - November 2023                                 | 50.9         | 37.9 | -0.71 |
| BIPEA Food Allergens – EGG 77B - Cake Mix – 09 - 0577 - November 2023                                 | 17.0         | 16.2 | -0.16 |
| BIPEA Food Allergens – EGG 77B - Cake Mix – 09 – 0677 - November 2023                                 | NOT DETECTED | <0.4 | -     |
| BIPEA Allergens in Industrial Process Water 77D - Industrial Process Water – 15 - 1077 – June 2023    | NOT DETECTED | <0.4 | -     |
| BIPEA Allergens in Industrial Process Water 77D - Industrial Process Water – 15 - 1277 – June 2023    | 23.1         | 22.3 | -0.09 |

|                                                                                                                                    |              |          |       |
|------------------------------------------------------------------------------------------------------------------------------------|--------------|----------|-------|
| BIPEA Food Allergens – EGG 77B - Pasta<br>– 08 - 2577 - March 2023                                                                 | 40.5         | 50.8     | 0.68  |
| BIPEA Food Allergens – EGG 77B - Pasta<br>– 07 - 2677 - March 2023                                                                 | NOT DETECTED | <0.4     | -     |
| BIPEA Food Allergens – EGG 77B - Pasta<br>– 07 - 2777 - March 2023                                                                 | 19.5         | 27.6     | 0.94  |
| BIPEA Allergens in Industrial Process<br>Water 77D - 14 2022-2023 - 0381<br>Industrial Process Water – 14 - 1077 -<br>January 2023 | 26.8         | 24.3     | -0.23 |
| BIPEA Allergens in Industrial Process<br>Water 77D - 14 2022-2023 - 0381<br>Industrial Process Water – 14 - 1177 -<br>January 2023 | 21.6         | 18.6     | -0.35 |
| BIPEA Allergens in Industrial Process<br>Water 77D - 14 2022-2023 - 0381<br>Industrial Process Water – 14 - 1277 -<br>January 2023 | NOT DETECTED | <0.4     | -     |
| BIPEA Food Allergens – EGG 77B - Cake<br>Mix – 08 – 0477 - October 2022                                                            | NOT DETECTED | <0.4     | -     |
| BIPEA Food Allergens – EGG 77B - Cake<br>Mix – 08 – 0577 - October 2022                                                            | 38.9         | 45.0     | 0.56  |
| BIPEA Food Allergens – EGG 77B - Cake<br>Mix – 08 – 0677 - October 2022                                                            | 11.6         | 19.8     | 1.67  |
| BIPEA Allergens in Industrial Process<br>Water 77D - Industrial Process Water –<br>13 - 1077 - June 2022                           | DETECTED     | DETECTED | -     |
| BIPEA Allergens in Industrial Process<br>Water 77D - Industrial Process Water –<br>13 - 1177 - June 2022                           | DETECTED     | DETECTED | -     |
| BIPEA Allergens in Industrial Process<br>Water 77D - Industrial Process Water –<br>13 - 1277 - June 2022                           | NOT DETECTED | <0.4     | -     |

| <b>Proficiency Test</b>                                                                               | <b>Assignment</b> | <b>Result (mg/kg)</b> | <b>Satisfactory results / results agreeing with assignment (%)</b> |
|-------------------------------------------------------------------------------------------------------|-------------------|-----------------------|--------------------------------------------------------------------|
| FAPAS Allergens on Cotton Swabs 27440A August – October 2025                                          | DETECTED          | DETECTED              | 97.0                                                               |
| FAPAS Allergens on Cotton Swabs 27440B August – October 2025                                          | NOT DETECTED      | NOT DETECTED          | 97.0                                                               |
| BIPEA Allergens in Industrial Process Water 77D - Industrial Process Water – 19 - 1077 – January 2025 | DETECTED          | DETECTED              | 100.0                                                              |
| BIPEA Allergens in Industrial Process Water 77D - Industrial Process Water – 18 - 1177 – January 2025 | DETECTED          | DETECTED              | 100.0                                                              |
| BIPEA Allergens in Industrial Process Water 77D - Industrial Process Water – 18 - 1277 – January 2025 | NOT DETECTED      | NOT DETECTED          | 100.0                                                              |
| BIPEA Allergens in Industrial Process Water 77D - Industrial Process Water – 17 - 1077 - June 2024    | DETECTED          | DETECTED              | 100.0                                                              |
| BIPEA Allergens in Industrial Process Water 77D - Industrial Process Water – 17 - 1177 - June 2024    | NOT DETECTED      | NOT DETECTED          | 100.0                                                              |
| BIPEA Allergens in Industrial Process Water 77D - Industrial Process Water – 17 - 1277 - June 2024    | DETECTED          | DETECTED              | 100.0                                                              |
| FAPAS Gluten, Egg and Milk on Swabs (Cotton) 27366E - September – November 2023                       | DETECTED          | DETECTED              | 97.0                                                               |

|                                                                                                    |              |              |       |
|----------------------------------------------------------------------------------------------------|--------------|--------------|-------|
| FAPAS Gluten, Egg and Milk on Swabs (Cotton) 27366F - September - November 2023                    | NOT DETECTED | NOT DETECTED | 97.0  |
| FAPAS Egg, Gluten and Milk on Swabs (Cotton) 27333A - September - November 2022                    | NOT DETECTED | NOT DETECTED | 93.0  |
| FAPAS Egg, Gluten and Milk on Swabs (Cotton) 27333B - September - November 2022                    | DETECTED     | DETECTED     | 97.0  |
| BIPEA Allergens in Industrial Process Water 77D - Industrial Process Water - 13 - 1077 - June 2022 | DETECTED     | DETECTED     | 100.0 |
| BIPEA Allergens in Industrial Process Water 77D - Industrial Process Water - 13 - 1177 - June 2022 | DETECTED     | DETECTED     | 100.0 |
| BIPEA Allergens in Industrial Process Water 77D - Industrial Process Water - 13 - 1277 - June 2022 | NOT DETECTED | <0.4         | 100.0 |

| <b>Proficiency Test</b>                                        | <b>Assigned value (mg/kg) according to R-Biopharm Ridascreen Fast Egg Ei/Egg protein R6402</b> | <b>Range for <math> z  \leq 2</math> according to R-Biopharm Ridascreen Fast Egg Ei/Egg protein R6402</b> | <b>Result (mg/kg)</b> | <b>Recovery (%)</b> |
|----------------------------------------------------------------|------------------------------------------------------------------------------------------------|-----------------------------------------------------------------------------------------------------------|-----------------------|---------------------|
| FAPAS Gluten, Milk & Egg in Cake Mix 27438A July - August 2025 | 18.7                                                                                           | 9.35 - 28.0                                                                                               | 15.8                  | 84.5                |
| FAPAS Allergens in Sausage Meat (Raw) 27427A March - May 2025  | 16.8                                                                                           | 8.42 - 25.3                                                                                               | 20.02                 | 119.2               |

|                                                                       |              |             |          |       |
|-----------------------------------------------------------------------|--------------|-------------|----------|-------|
| FAPAS Allergens in Sausage Meat (Raw) 27427B March – May 2025         | NOT DETECTED | -           | <1.0     | -     |
| FAPAS Gluten, Milk & Egg in Cooked Biscuit 27418 January – March 2025 | 5.33         | 2.67 – 8.00 | 3.45     | 64.7  |
| FAPAS Allergens in Cooked Biscuit 27382 January – March 2024          | 61.5         | 30.7 – 92.2 | 47.6     | 77.4  |
| FAPAS Allergens in Cake Mix - 27377A - November 2023 – January 2024   | 29.3         | 14.6 – 43.9 | 33.4     | 114.0 |
| FAPAS Allergens in Cake Mix - 27377B - November 2023 – January 2024   | 36.1         | 18.0 – 54.1 | 38.6     | 106.9 |
| FAPAS Allergens in Cake Mix - 27364A - July – September 2023          | 17.6         | 8.81 – 26.4 | 14.1     | 80.1  |
| FAPAS Allergens in Cake Mix - 27364B - July – September 2023          | 33.6         | 16.8 – 50.4 | 28.7     | 85.4  |
| FAPAS Allergens in Cake Mix - 27331A - July – September 2022          | 5.66         | 2.83 – 8.48 | 5.40     | 95.4  |
| FAPAS Allergens in Cake Mix - 27331B - July – September 2022          | 7.16         | 3.58 – 10.7 | 7.41     | 103.5 |
| FAPAS Allergens in Cooked Biscuit - 27314 - January – March 2022      | 21.8         | 10.9 – 32.6 | 18.3     | 83.9  |
| FAPAS Allergens in Cake Mix - 27311A - November 2021 – February 2022  | DETECTED     | -           | DETECTED | -     |

|                                                                            |              |              |          |       |
|----------------------------------------------------------------------------|--------------|--------------|----------|-------|
| FAPAS Allergens in Cake Mix<br>- 27311B - November 2021 –<br>February 2022 | DETECTED     | -            | DETECTED | -     |
| FAPAS Allergens in Cooked<br>Biscuit – 27286 - January –<br>March 2021     | 13.0         | 6.5 – 19.5   | 15.0     | 115.4 |
| FAPAS Egg in Cake Mix -<br>27296A - April – June 2021                      | NOT DETECTED | -            | <0.2     | -     |
| FAPAS Egg in Cake Mix -<br>27296B - April – June 2021                      | 15.38        | 6.43 – 19.3  | 12.9     | 119.2 |
| FAPAS Allergens in Cake Mix<br>- 27301A - July – September<br>2021         | 9.03         | 4.51 – 13.50 | 5.34     | 59.1  |
| FAPAS Allergens in Cake Mix<br>- 27301B - July – September<br>2021         | 6.72         | 3.36 – 10.10 | 3.55     | 52.9  |

| <b>Proficiency Test</b>                                                                                        | <b>Assignment</b> | <b>Result (mg/kg)</b> | <b>Satisfactory results /<br/>results agreeing with<br/>assignment (%)</b> |
|----------------------------------------------------------------------------------------------------------------|-------------------|-----------------------|----------------------------------------------------------------------------|
| BIPEA Allergens in Industrial<br>Process Water 77D - Industrial<br>Process Water – 19 - 1077 – January<br>2025 | DETECTED          | DETECTED              | 100.0                                                                      |
| BIPEA Allergens in Industrial<br>Process Water 77D - Industrial<br>Process Water – 18 - 1177 – January<br>2025 | DETECTED          | DETECTED              | 100.0                                                                      |
| BIPEA Allergens in Industrial<br>Process Water 77D - Industrial<br>Process Water – 18 - 1277 – January<br>2025 | NOT DETECTED      | NOT DETECTED          | 100.0                                                                      |
| BIPEA Allergens in Industrial<br>Process Water 77D - Industrial                                                | DETECTED          | DETECTED              | 100.0                                                                      |

|                                                                                                    |              |              |       |
|----------------------------------------------------------------------------------------------------|--------------|--------------|-------|
| Process Water – 17 - 1077 - June 2024                                                              |              |              |       |
| BIPEA Allergens in Industrial Process Water 77D - Industrial Process Water – 17 - 1177 - June 2024 | NOT DETECTED | NOT DETECTED | 100.0 |
| BIPEA Allergens in Industrial Process Water 77D - Industrial Process Water – 17 - 1277 - June 2024 | DETECTED     | DETECTED     | 100.0 |
| FAPAS Gluten, Egg and Milk on Swabs (Cotton) 27366E - September – November 2023                    | DETECTED     | DETECTED     | 97.0  |
| FAPAS Gluten, Egg and Milk on Swabs (Cotton) 27366F - September – November 2023                    | NOT DETECTED | NOT DETECTED | 97.0  |
| FAPAS Egg, Gluten and Milk on Swabs (Cotton) 27333A - September – November 2022                    | NOT DETECTED | NOT DETECTED | 93.0  |
| FAPAS Egg, Gluten and Milk on Swabs (Cotton) 27333B - September – November 2022                    | DETECTED     | DETECTED     | 97.0  |
| BIPEA Allergens in Industrial Process Water 77D - Industrial Process Water – 13 - 1077 - June 2022 | DETECTED     | DETECTED     | 100.0 |
| BIPEA Allergens in Industrial Process Water 77D - Industrial Process Water – 13 - 1177 - June 2022 | DETECTED     | DETECTED     | 100.0 |
| BIPEA Allergens in Industrial Process Water 77D - Industrial Process Water – 13 - 1277 - June 2022 | NOT DETECTED | <0.4         | 100.0 |

