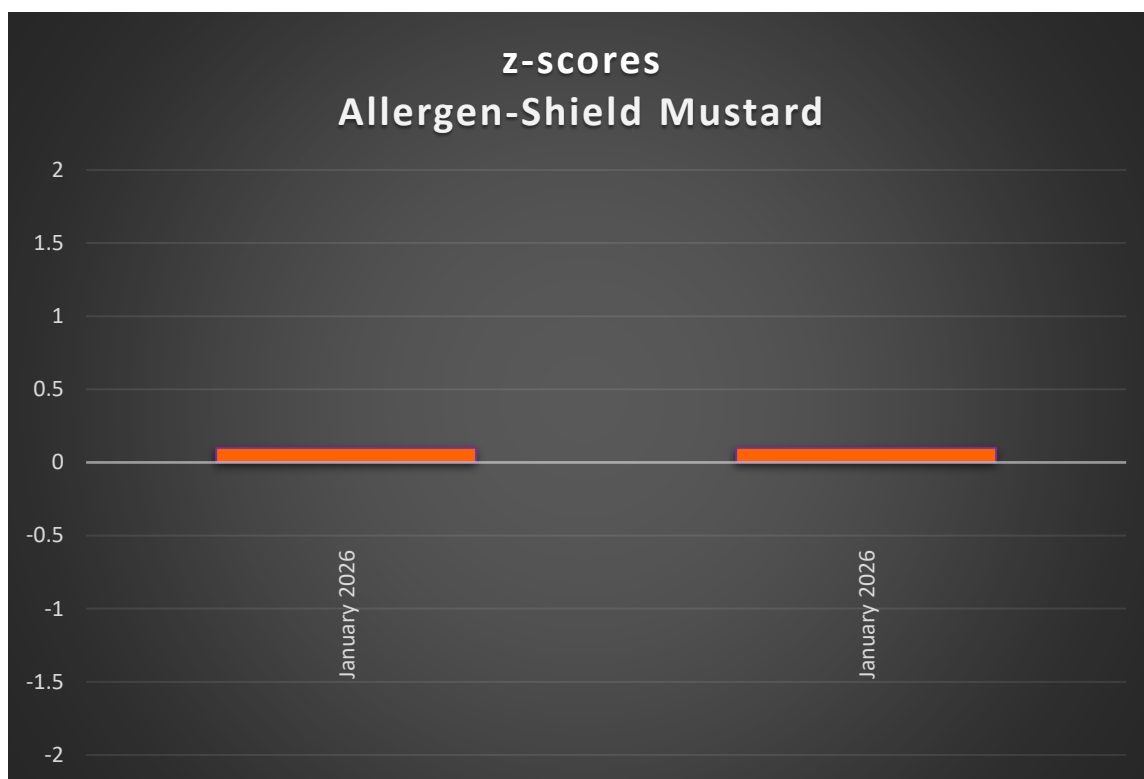




Proficiency Test	Assigned value (mg/L)	Result (mg/L)	Z-score
FAPAS - Mustard in Durum Wheat – 27456A – December 2025 – January 2026	51.5	52.9	0.1
FAPAS - Mustard in Durum Wheat – 27456B – December 2025 – January 2026	104	106	0.1

Proficiency Test	Assignment	Result (mg/kg)	Satisfactory results / results agreeing with assignment (%)
FAPAS Celery and Mustard in Instant Soup Powder – 27474A – June – July 2026	DETECTED	DETECTED	98.0

FAPAS Celery and Mustard in Instant Soup Powder – 27474B - June – July 2026	DETECTED	DETECTED	93.0
FAPAS Celery and Mustard in Instant Soup Powder – 27434A - June – July 2025	DETECTED	DETECTED	75.0
FAPAS Celery and Mustard in Instant Soup Powder – 27434B - June – July 2025	DETECTED	DETECTED	100.0

Proficiency Test	Assigned value (mg/kg) according to R-Biopharm Ridascreen Fast Senf/Mustard R6152	Range for $z \leq 2$ according to R-Biopharm Ridascreen Fast Senf/Mustard R6152	Result (mg/kg)	Recovery (%)
FAPAS Celery & Mustard in Instant Soup Powder 27397A June – July 2024	43.4	21.7 – 65.1	41.5	95.6
FAPAS Celery & Mustard in Instant Soup Powder 27397B June – July 2024	NOT DETECTED	-	<1.0	-
FAPAS Celery & Mustard in Instant Soup Powder 27361A June – July 2023	44.1	22.0 – 66.1	46.1	104.5
FAPAS Celery & Mustard in Instant Soup Powder 27361B June – July 2023	NOT DETECTED	-	<1.0	-
FAPAS Celery & Mustard in Instant Soup Powder 27328A June – August 2022	NOT DETECTED	-	<1.0	-
FAPAS Celery & Mustard in Instant Soup Powder 27328B June – August 2022	DETECTED	-	DETECTED	-

