




Proficiency Test	Assigned value (mg/L)	Result (mg/L)	Z-score
FAPAS - Almond In Cake - 27436A - June – August 2025 (Allergenic Protein)	4.41	4.41	0.0
FAPAS - Almond In Cake – 27436B - June – August 2025 (Allergenic Protein)	NOT DETECTED	<0.21	-
BIPEA - Food Allergens - Nuts - 77G - Yogurt – 1 – 4677 - December 2024	58.5	52.8	-0.67
BIPEA - Food Allergens - Nuts - 77G - Yogurt – 1 – 4577 - December 2024	NOT DETECTED	<2.0	-
BIPEA - Food Allergens - Nuts - 77G - Yogurt – 1 – 4477 - December 2024	112.4	80.5	-1.48

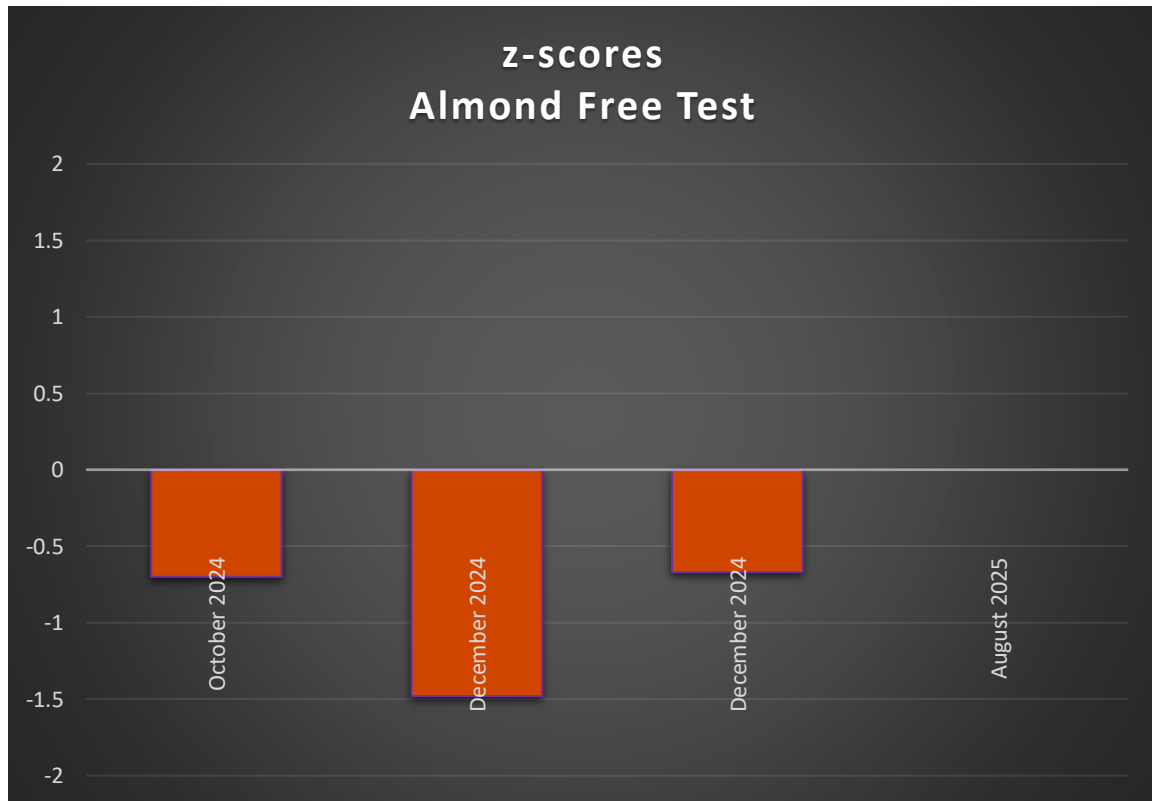

FAPAS - Almond In Chocolate - 27406A - September – October 2024 (Allergenic Protein)	4.28	3.57	-0.7
FAPAS - Almond In Chocolate – 27406B - September – October 2024	NOT DETECTED	<2.0	-

Proficiency Test	Assigned value (mg/kg) according to Allergen-Shield Almond A1548/A1596	Range for $ z \leq 2$ according to Allergen-Shield Almond A1548/A1596	Result (mg/kg)	Recovery (%)
FAPAS - Almond In Chocolate - 27406A - September – October 2024	17.8	8.88 – 26.6	17.0	95.5

Proficiency Test	Assignment	Result (mg/kg)	Satisfactory results / results agreeing with assignment (%)
FAPAS - Nuts In Chocolate - 27414 - November 2024 – January 2025	NOT DETECTED	NOT DETECTED	95.0
FAPAS - Nuts In Chocolate - 27378 - November 2023 – January 2024	DETECTED	DETECTED	100.0