




Proficiency Test	Assigned value (mg/L)	Result (mg/L)	Z-score
BIPEA Food Allergens - Milk 77C - 20 - Infant Formula - 21 - 0777 - April 2025	NOT DETECTED	<0.1	-
BIPEA Food Allergens - Milk 77C - 20 - Infant Formula - 21 - 0877 - April 2025	1.2	0.8	-1.11
BIPEA Food Allergens - Milk 77C - 20 - Infant Formula - 20 - 0977 - April 2025	0.5	0.4	-0.50
BIPEA Food Allergens - Milk 77C - 19 - Infant Formula - 20 - 0777 - November 2024	NOT DETECTED	<0.6	-

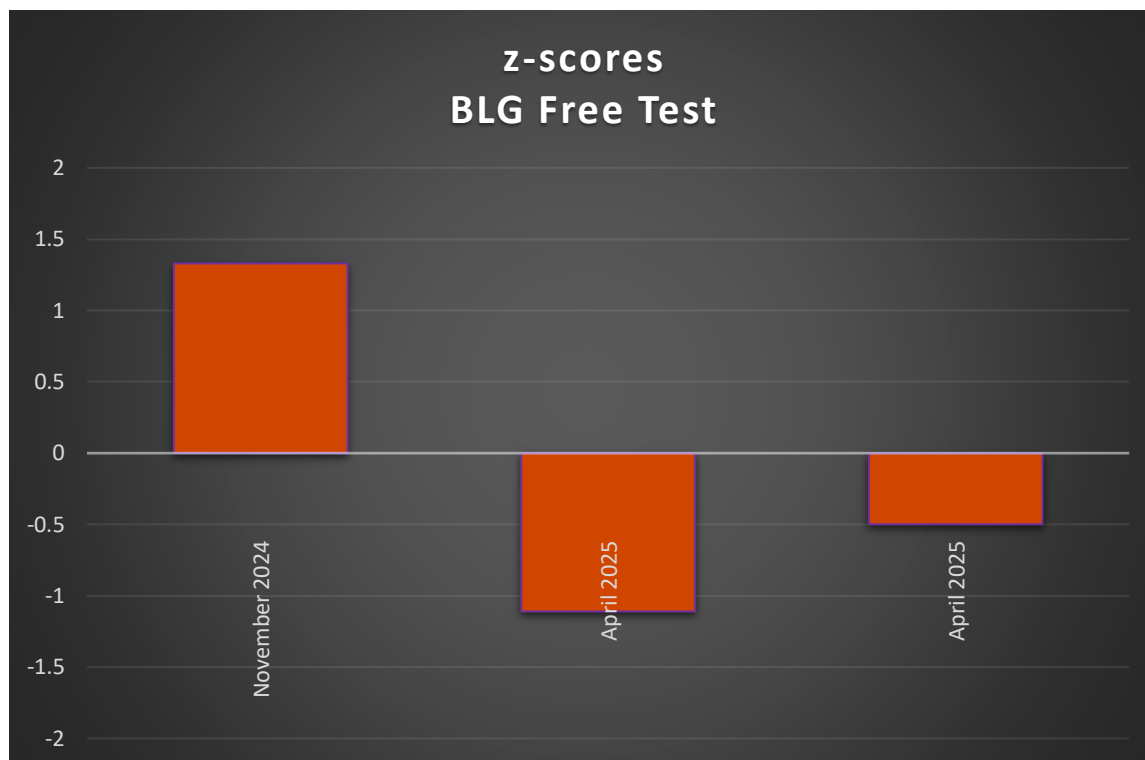

BIPEA Food Allergens - Milk 77C - 19 - Infant Formula - 20 - 0877 - November 2024	NO VALUE WAS ASSIGNED	<0.6	-
BIPEA Food Allergens - Milk 77C - 19 - Infant Formula - 19 - 0977 - November 2024	0.4	0.6	1.33

Proficiency Test	Assigned value (mg/kg) according to R-Biopharm Ridascreen Fast β -Lactoglobulin R4912	Range for $ z \leq 2$ according to R-Biopharm Ridascreen Fast Fast β -Lactoglobulin R4912	Result (mg/kg)	Recovery (%)
FAPAS Milk in Infant Breakfast Cereal 27435A June - August 2025	1.72	0.858 - 2.57	2.26	131.4
FAPAS Milk in Infant Breakfast Cereal 27435B June - August 2025	NOT DETECTED	-	<0.1	-

Proficiency Test	Assignment	Result (mg/kg)	Satisfactory results / results agreeing with assignment (%)
FAPAS Milk & Gluten in Infant Soya Formula 27422A February - March 2025	DETECTED	DETECTED	-
FAPAS Milk & Gluten in Infant Soya Formula 27422B February - March 2025	NOT DETECTED	<1.50	-
FAPAS Gluten, Milk & Egg in Cooked Biscuit 27418 January - March 2025	DETECTED	DETECTED	-



FAPAS Milk in Infant Breakfast Cereal - 27398A - June - August 2024	NOT DETECTED	<2.0	-
FAPAS Milk in Infant Breakfast Cereal - 27398B - June - August 2024	DETECTED	DETECTED	-