




Proficiency Test	Assigned value (mg/L)	Result (mg/L)	Z-score
BIPEA Food Allergens - Milk 77C - 20 - Infant Formula - 21 - 0777 - April 2025	NOT DETECTED	<5.0	-
BIPEA Food Allergens - Milk 77C - 20 - Infant Formula - 21 - 0877 - April 2025	26.6	40.0	1.44
BIPEA Food Allergens - Milk 77C - 20 - Infant Formula - 20 - 0977 - April 2025	10.1	11.0	0.43
BIPEA Food Allergens - Milk 77C - 19 - Infant Formula - 20 - 0777 - November 2024	NOT DETECTED	<4.5	-

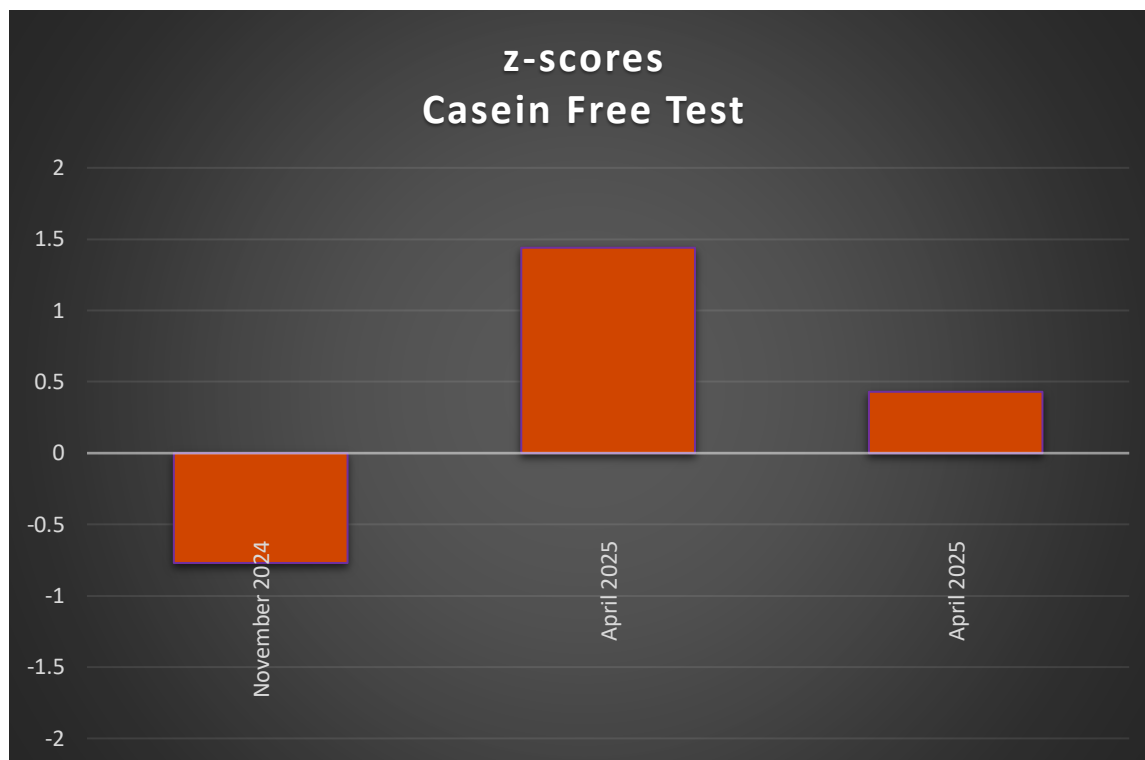

BIPEA Food Allergens - Milk 77C - 19 - Infant Formula - 20 - 0877 - November 2024	2.0	<4.5	-
BIPEA Food Allergens - Milk 77C - 19 - Infant Formula - 19 - 0977 - November 2024	15.7	12.5	-0.77

Proficiency Test	Assigned value (mg/kg) according to R-Biopharm Ridascreen Fast Casein R4612	Range for $ z \leq 2$ according to R-Biopharm Ridascreen Fast Casein R4612	Result (mg/kg)	Recovery (%)
FAPAS Milk in Infant Breakfast Cereal 27435A June - August 2025	14.0	7.01 - 21.0	15.4	110.0
FAPAS Milk in Infant Breakfast Cereal 27435B June - August 2025	NOT DETECTED	-	<5.0	-
FAPAS Milk & Gluten in Infant Soya Formula 27422A February - March 2025	19.7	9.86 - 29.6	16.0	81.2
FAPAS Milk & Gluten in Infant Soya Formula 27422B February - March 2025	NOT DETECTED	-	<6.0	-
FAPAS Gluten, Milk & Egg in Cooked Biscuit 27418 January - March 2025	17.9	8.93 - 26.8	14.6	81.6

Proficiency Test	Assigned value (mg/kg) according to Allergen-Shield	Range for $ z \leq 2$ according to Allergen-Shield Casein A1248/A1296	Result (mg/kg)	Recovery (%)
------------------	---	--	----------------	--------------



	Casein A1248/A1296			
FAPAS Gluten, Egg and Milk in Cake Mix 27400A - July - September 2024	14.4	7.18 - 21.5	14.0	97.2
FAPAS Gluten, Egg and Milk in Cake Mix 27400B - July - September 2024	30.4	15.2 - 45.5	30.0	98.7

Proficiency Test	Assigned value (mg/kg) according to R-Biopharm Ridascreen Fast Casein R4612	Range for $ z \leq 2$ according to R-Biopharm Ridascreen Fast Casein R4612	Result (mg/kg)	Recovery (%)
FAPAS Milk in Infant Breakfast Cereal - 27398A - June - August 2024	NOT DETECTED	-	<4.0	-
FAPAS Milk in Infant Breakfast Cereal - 27398B - June - August 2024	19.1	9.55 - 28.6	12.8	67.0