




Proficiency Test	Assigned value (mg/L)	Result (mg/L)	Z-score
FAPAS Milk in Infant Breakfast Cereal 27475A June – August 2026	18.1		
FAPAS Milk in Infant Breakfast Cereal 27475B June – August 2026	NOT DETECTED	<5.0	-
BIPEA Allergens in Industrial Process Water 77D - Industrial Process Water – 21 - 1077 – June 2026	17.8	11.5	-0.91
BIPEA Allergens in Industrial Process Water 77D - Industrial Process Water – 21 - 1177 – June 2026	NOT DETECTED	<5.0	-
BIPEA Allergens in Industrial Process Water 77D - Industrial Process Water – 21 - 1277 – June 2026	24.8	16.2	-1.18
BIPEA Food Allergens - Milk 77C - Infant Formula – 23 - 0777 – April 2026	6.7	8.2	0.64

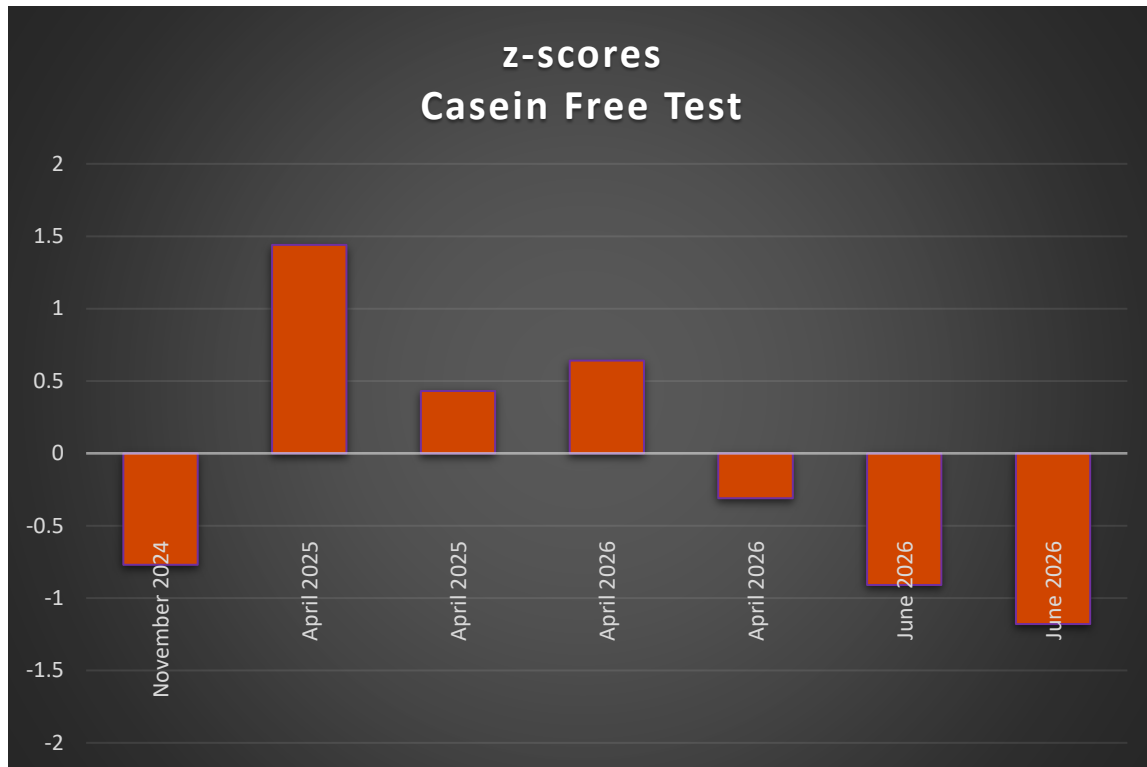

BIPEA Food Allergens - Milk 77C - Infant Formula – 23 - 0877 – April 2026	NOT DETECTED	<5.0	-
BIPEA Food Allergens - Milk 77C - Infant Formula – 22 - 0977 – April 2026	20.2	18.0	-0.31
BIPEA Food Allergens - Milk 77C - 20 - Infant Formula – 21 - 0777 – April 2025	NOT DETECTED	<5.0	-
BIPEA Food Allergens - Milk 77C - 20 - Infant Formula – 21 - 0877 – April 2025	26.6	40.0	1.44
BIPEA Food Allergens - Milk 77C - 20 - Infant Formula – 20 - 0977 – April 2025	10.1	11.0	0.43
BIPEA Food Allergens - Milk 77C - 19 - Infant Formula – 20 - 0777 – November 2024	NOT DETECTED	<4.5	-
BIPEA Food Allergens - Milk 77C - 19 - Infant Formula – 20 - 0877 – November 2024	2.0	<4.5	-
BIPEA Food Allergens - Milk 77C - 19 - Infant Formula – 19 - 0977 – November 2024	15.7	12.5	-0.77

Proficiency Test	Assigned value (mg/kg) according to Allergen-Shield Casein A1248/A1296	Range for $ z \leq 2$ according to Allergen-Shield Casein A1248/A1296	Result (mg/kg)	Recovery (%)
FAPAS Milk in Infant Breakfast Cereal 27475A June – August 2026	18.1	9.05 – 27.2	17.3	95.6
FAPAS Milk in Infant Breakfast Cereal 27475B June – August 2026	NOT DETECTED	-	<5.0	-
FAPAS Milk & Gluten in Infant Soya Formula	NOT DETECTED	-	<5.0	-

27463A March - May 2026				
FAPAS Milk & Gluten in Infant Soya Formula 27463B March - May 2026	18.1	9.05 - 27.2	13.68	75.6
FAPAS Gluten, Milk & Egg in Cooked Biscuit 27459 January - March 2026	19.9	9.94 - 29.8	25.96	130.5
FAPAS Gluten, Egg and Milk in Cake Mix 27400A - July - September 2024	14.4	7.18 - 21.5	14.0	97.2
FAPAS Gluten, Egg and Milk in Cake Mix 27400B - July - September 2024	30.4	15.2 - 45.5	30.0	98.7

Proficiency Test	Assigned value (mg/kg) according to R-Biopharm Ridascreen Fast Casein R4612	Range for $z \leq 2$ according to R-Biopharm Ridascreen Fast Casein R4612	Result (mg/kg)	Recovery (%)
FAPAS Gluten, Milk & Egg in Cake Mix 27438A July - August 2025	17.5	8.74 - 26.2	17.0	97.1
FAPAS Gluten, Milk & Egg in Cake Mix 27438B July - August 2025	26.5	13.3 - 39.8	31.5	118.9
FAPAS Milk in Infant Breakfast Cereal 27435A June - August 2025	14.0	7.01 - 21.0	15.4	110.0
FAPAS Milk in Infant Breakfast Cereal 27435B June - August 2025	NOT DETECTED	-	<5.0	-



FAPAS Milk & Gluten in Infant Soya Formula 27422A February – March 2025	19.7	9.86 – 29.6	16.0	81.2
FAPAS Milk & Gluten in Infant Soya Formula 27422B February – March 2025	NOT DETECTED	-	<6.0	-
FAPAS Gluten, Milk & Egg in Cooked Biscuit 27418 January – March 2025	17.9	8.93 – 26.8	14.6	81.6

Proficiency Test	Assigned value (mg/kg) according to R- Biopharm Ridascreen Fast Casein R4612	Range for $ z \leq 2$ according to R- Biopharm Ridascreen Fast Casein R4612	Result (mg/kg)	Recovery (%)
FAPAS Milk in Infant Breakfast Cereal - 27398A - June - August 2024	NOT DETECTED	-	<4.0	-
FAPAS Milk in Infant Breakfast Cereal - 27398B - June - August 2024	19.1	9.55 – 28.6	12.8	67.0