




Proficiency Test	Assigned value (mg/L or mg/kg)	Result (mg/L or mg/kg)	Z-score
FAPAS Gluten, Milk & Egg in Cooked Biscuit 27418 January – March 2025	22.2	16.5	-1.0
BIPEA Food Allergens 77A - Gluten – Cake 7 - 3477 – February 2025	NOT DETECTED	<5.0	-
BIPEA Food Allergens 77A - Gluten – Cake 7 - 3577 – February 2025	32.0	20.0	-1.33
BIPEA Food Allergens 77A - Gluten – Cake 7 - 3677 – February 2025	66.0	40.0	-1.63

Proficiency Test	Assigned value (mg/kg) according to Allergen-Shield Gluten S A1048/A1096	Range for $ z \leq 2$ according to Allergen-Shield Gluten S A1048/A1096	Result (mg/kg)	Recovery (%)
FAPAS Gluten, Egg and Milk in Cake Mix 27400A July – September 2024	12.8	6.38 – 19.1	12.0	93.8
FAPAS Gluten, Egg and Milk in Cake Mix 27400B July – September 2024	24.3	12.1 – 36.4	22.0	90.5
BIPEA Food Allergens 77A - Gluten – Infant Flour 13 - 0177 - September 2024	24.0	8.00 – 40.0	30.0	125.0
BIPEA Food Allergens 77A - Gluten – Infant Flour 11 - 0277 - September 2024	79.0	51.0 – 107	80.0	101.3
BIPEA Food Allergens 77A - Gluten – Infant Flour 12 - 0377 - September 2024	NOT DETECTED	-	<5.0	-

Proficiency Test	Assigned value (mg/kg) according to R-Biopharm Ridascreen Gliadin R7001	Range for $ z \leq 2$ according to R- Biopharm Ridascreen Gliadin R7001	Result (mg/kg)	Recovery (%)
FAPAS Milk & Gluten in Infant Soya Formula 27422A February – March 2025	37.5	18.7 – 56.2	28.6	76.3
FAPAS Milk & Gluten in Infant Soya Formula 27422B February – March 2025	NOT DETECTED	-	<5.0	-
FAPAS Allergens in Sausage Meat (Raw) 27427A March – May 2025	18.5	9.23 – 27.7	20.2	109.2

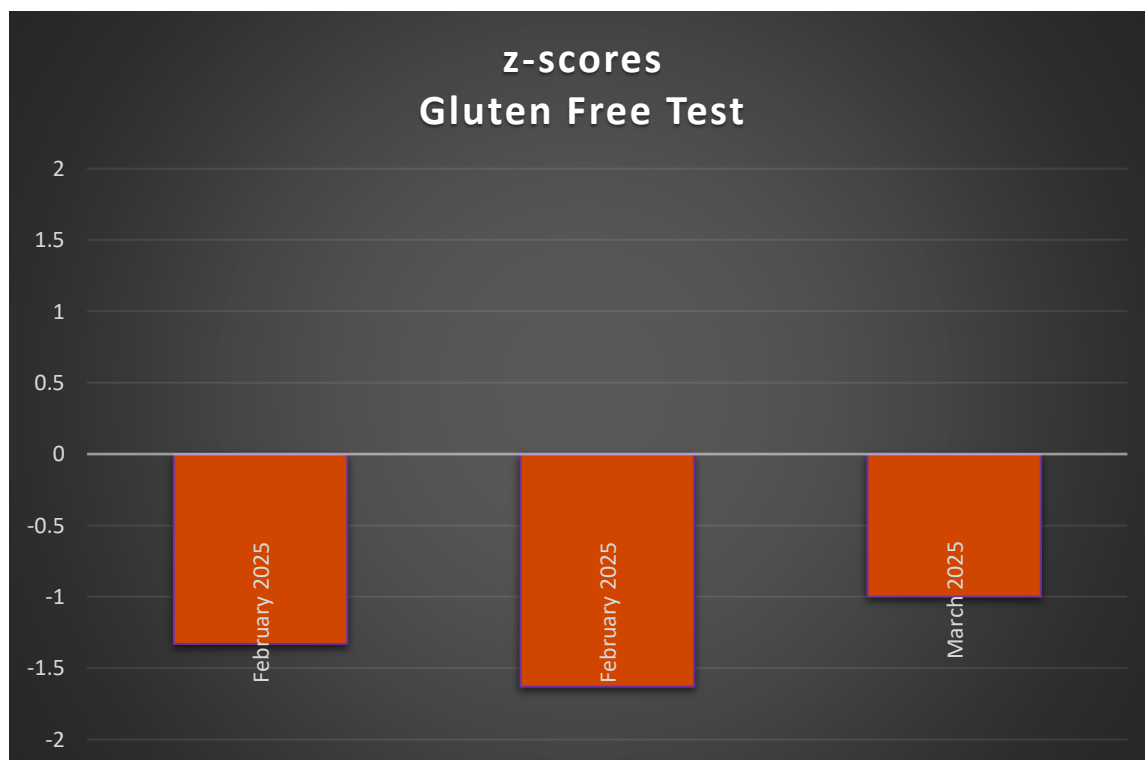

FAPAS Allergens in Sausage Meat (Raw) 27427B March – May 2025	NOT DETECTED	-	<5.0	-
FAPAS Gluten & Sesame in Cumin 27423A February – April 2025	50.1	25.1 – 75.2	52.5	104.8
FAPAS Gluten & Sesame in Cumin 27423B February – April 2025	NOT DETECTED	-	<3.0	-

Proficiency Test	Assigned value (mg/kg)	Range for $ z \leq 2$	Result (mg/kg)	Recovery (%)
BIPEA Food Allergens 77A - Gluten – Infant Flour 13 - 0177 - September 2024	24.0	8.00 – 40.0	30.0	125.0
BIPEA Food Allergens 77A - Gluten – Infant Flour 11 - 0277 - September 2024	79.0	51.0 – 107	80.0	101.3
BIPEA Food Allergens 77A - Gluten – Infant Flour 12 - 0377 - September 2024	NOT DETECTED	-	<5.0	-