




Proficiency Test	Assigned value (mg/L)	Result (mg/L)	Z-score
BIPEA Food Allergens - Soya - 77H - Wheat Flour - 09 - 3177 - June 2025	50.4	60.6	1.07
BIPEA Food Allergens - Soya - 77H - Wheat Flour - 09 - 3277 - June 2025	NOT DETECTED	<1.5	-
BIPEA Food Allergens - Soya - 77H - Wheat Flour - 10 - 3377 - June 2025	18.3	18.8	0.33
FAPAS Soya in Wheat Flour 27392B - April - May 2024	21.9	25.0	0.6

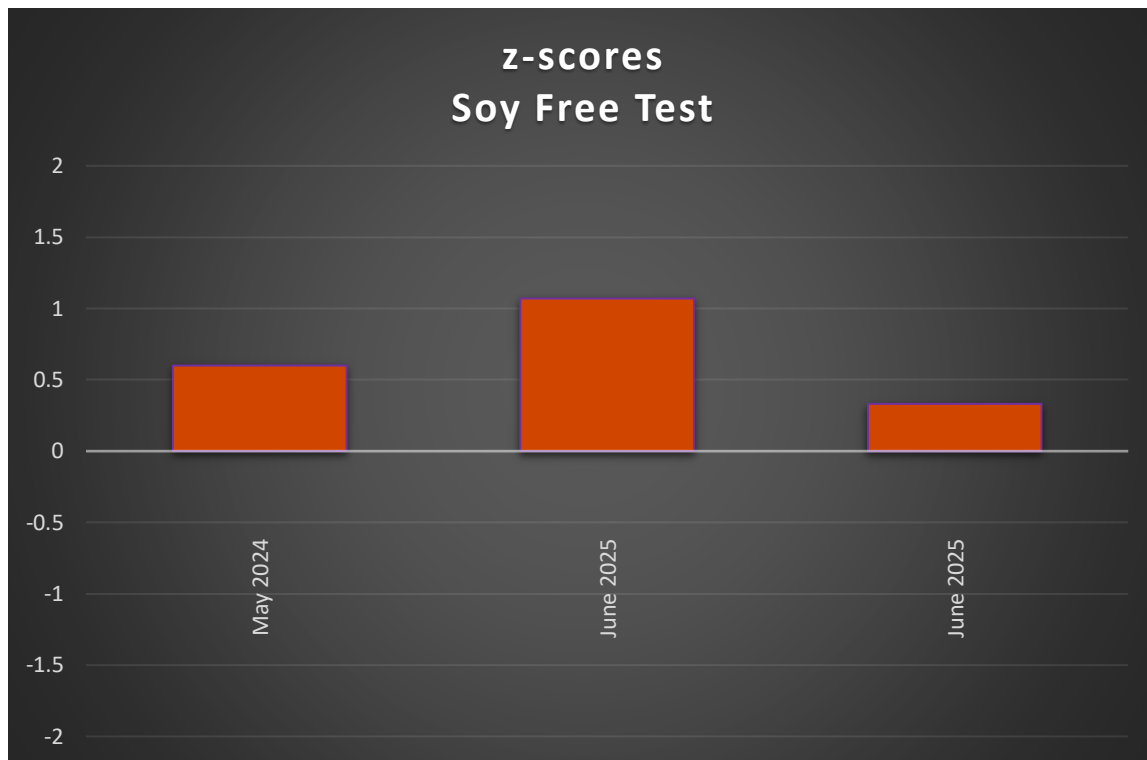

Proficiency Test	Assignment	Result (mg/kg)	Satisfactory results / results agreeing with assignment (%)
FAPAS Allergens on Swab 27419E - May - July 2025	NOT DETECTED	<0.1	91.0
FAPAS Allergens on Swab 27419F - May - July 2025	DETECTED	DETECTED	39.0
BIPEA Food Allergens - Soya - 77H - Meat - 10 - 2277 - January 2025	DETECTED	DETECTED	92.0
BIPEA Food Allergens - Soya - 77H - Meat - 10 - 2377 - January 2025	NOT DETECTED	NOT DETECTED	92.0
BIPEA Food Allergens - Soya - 77H - Meat - 10 - 2477 - January 2025	DETECTED	DETECTED	92.0
FAPAS Soya in Wheat Flour 27392A - April - May 2024	DETECTED	DETECTED	100.0
FAPAS Soya in Wheat Flour 27392B - April - May 2024	DETECTED	DETECTED	100.0
FAPAS Allergens on Swab 27383E - January - March 2024	NOT DETECTED	NOT DETECTED	100.0
FAPAS Allergens on Swab 27383F - January - March 2024	DETECTED	DETECTED	63.0